

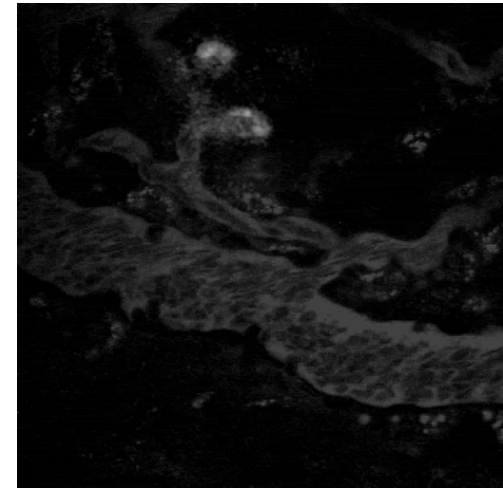


# Building Your Own 2-Photon Microscope: Challenges, Advantages and Limitations

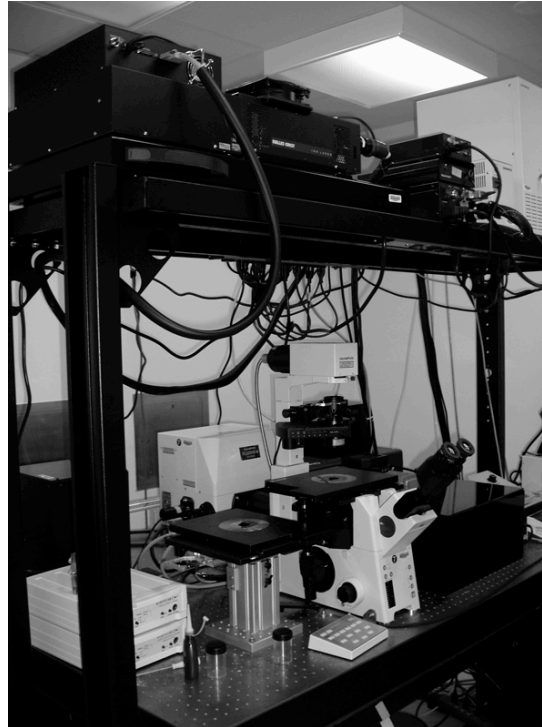

**Roberto Weigert, Ph.D.**  
**Intracellular Membrane Trafficking Unit**  
**Oral and Pharyngeal Cancer Branch**  
**NIDCR-NIH**

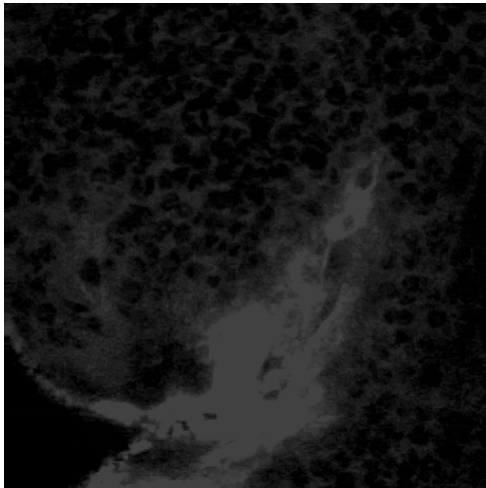

# Building Your Own 2-Photon Microscope: Challenges, Advantages and Limitations

How did we manage to build a 2-photon microscope?

12% Ethanol



A phone:

317-278-0436

Internet connection:

[kwdunn@iupui.edu](mailto:kwdunn@iupui.edu)

Espresso Machine  
“Gaggia”



# Building Your Own 2-Photon Microscope

Turn Key System

Why?

Build your own

More expensive  
\$400K-500K

\$150K Laser

- 1) Budget
- 1) Initial expenses
- 2) Maintenance
- 2) Flexibility
- 3) Upgrades

Buy Confocal microscope



Convert to a 2-photon

First step (first major decision)

Which Platform?

Olympus

- 1) We copied the system built here
- 2) Flexibility
- 3) Support

# Upright vs. Inverted

Flexibility:

- 1) 2-photon
- 2) Confocal

- 1) Intravital imaging
- 2) Live Cell imaging




- 1) All the organs

## Motion artifacts

**“Custom made holding device specifically designed for the organ of interest”**

w/o holding device

with holding device

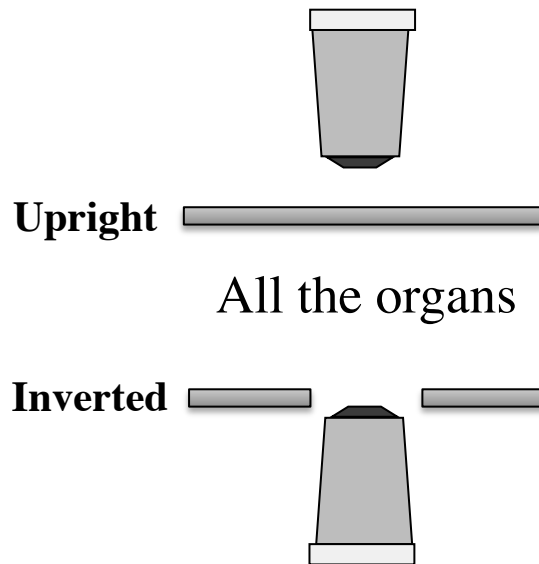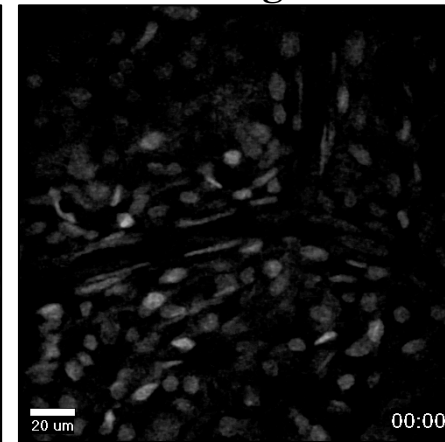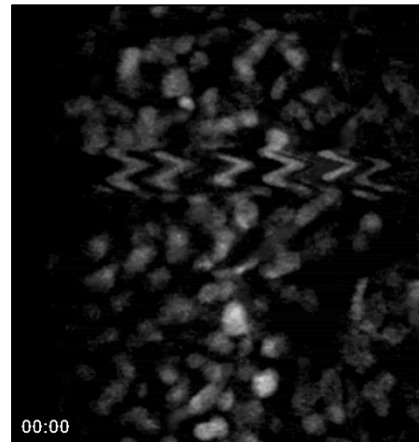

**Upright**




All the organs

**Inverted**


All the organs but the brain  
Live Cell Imaging

**“Positioning and securing the organ to the coverslip”**


# Upright vs. Inverted

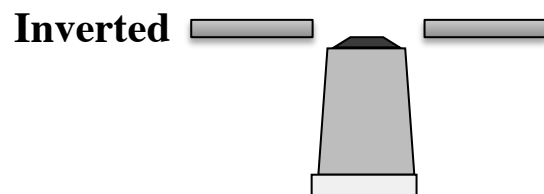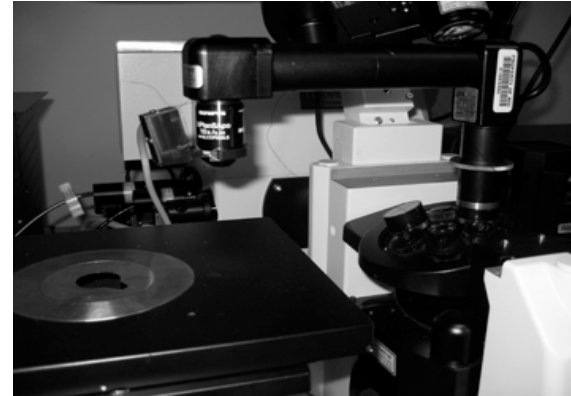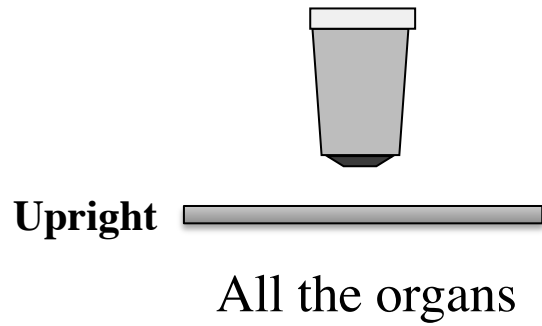

**Inverted  
converted to  
upright**

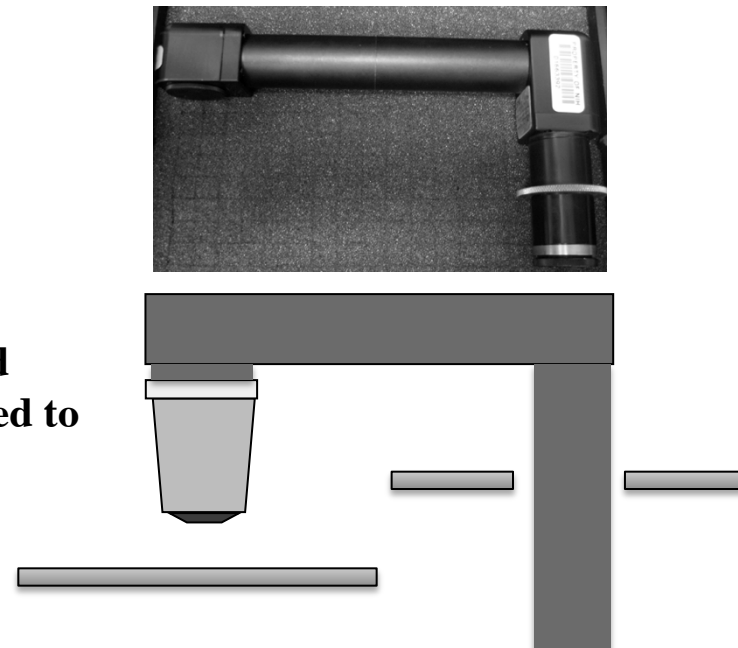

# Upright vs. Inverted

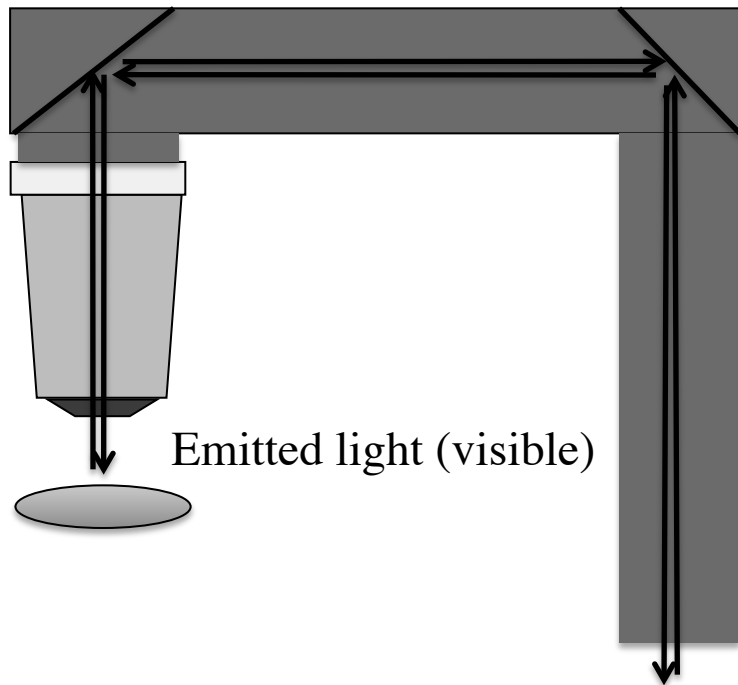


Excitation beam (IR)

- 1) Optimized for visible light
- 2) Increase the light path
  - 1) Model available with PMT on top
- 3) Loss of power (5-10%)
- 4) No effects on laser pulse width
- 5) Requires extra stage
- 6) Head can be rotated
- 7) Adaptors for lenses



# Upright vs. Inverted

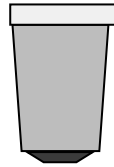


**Upright**




All the organs

**Upright  
converted to  
inverted**

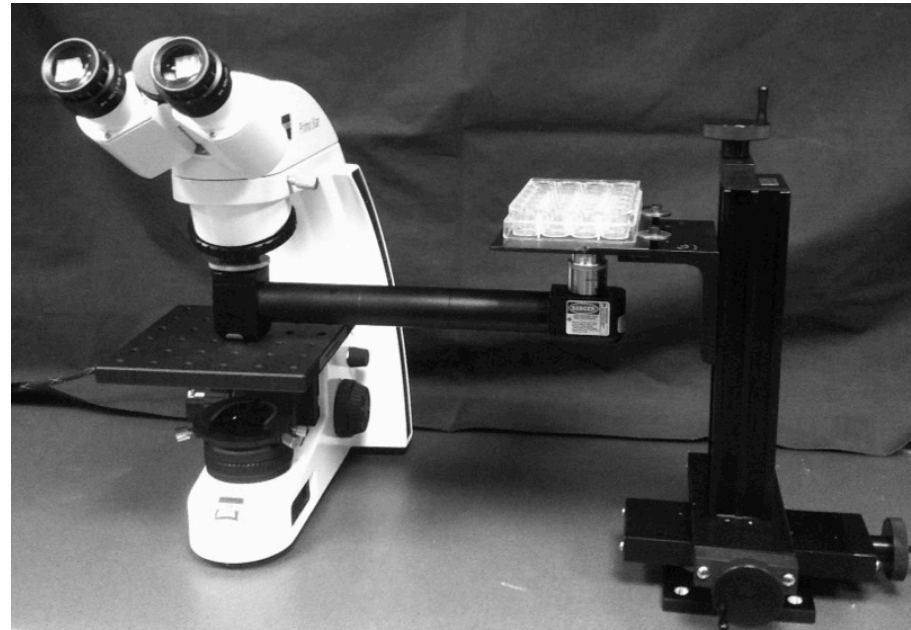

# Laser

**488 nm**  
**Laser combiner 561 nm**  
**633 nm**  
**UV laser 405 nm**

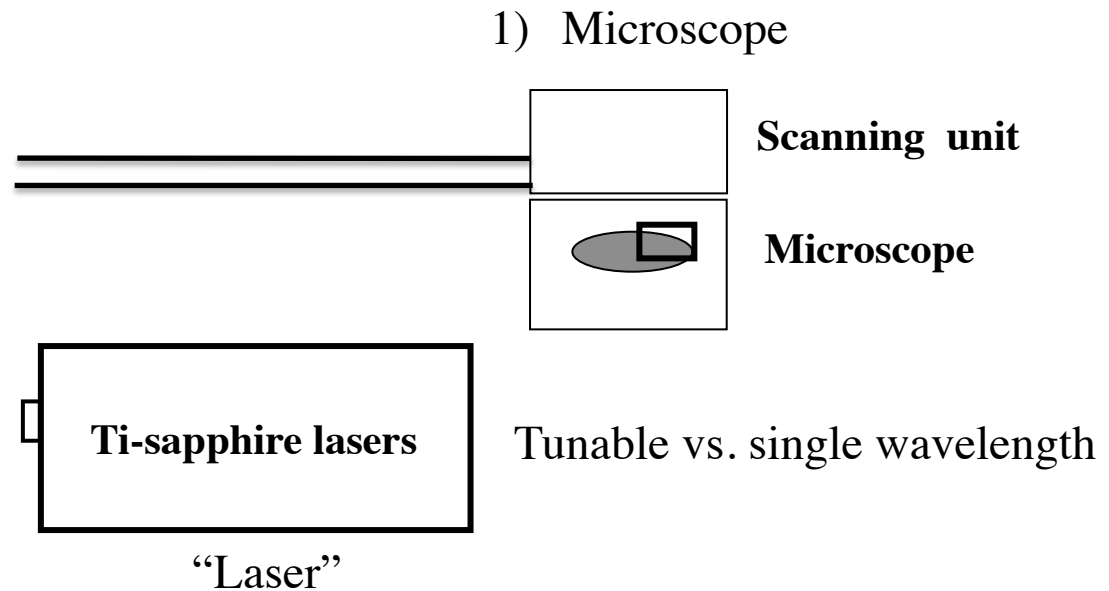

**High power lasers (3-4 W)**  
**Repetition rates: 80-100 Hz**


**Pulses: 100-150 fs**

**Beam diameter: 1.2 +/- 0.2 mm**




**Tunable: 680-1080 nm**




# Laser output power

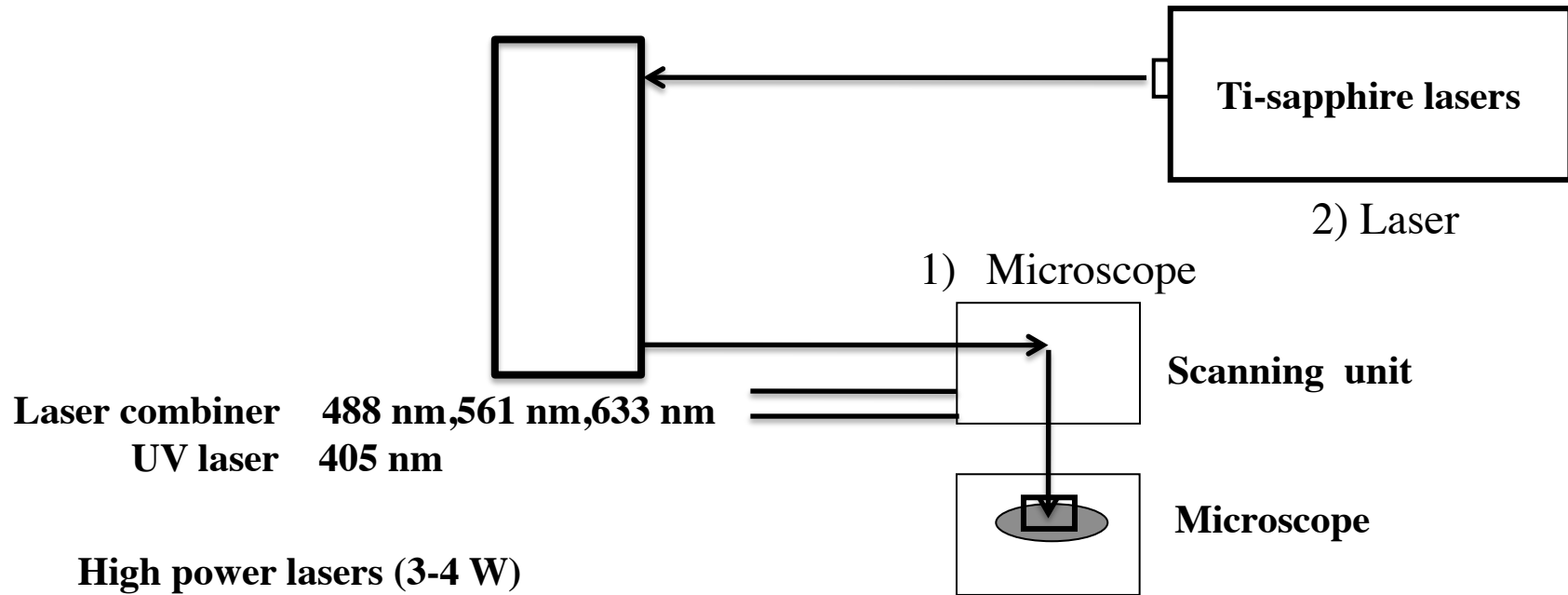


## High power lasers (3-4 W)

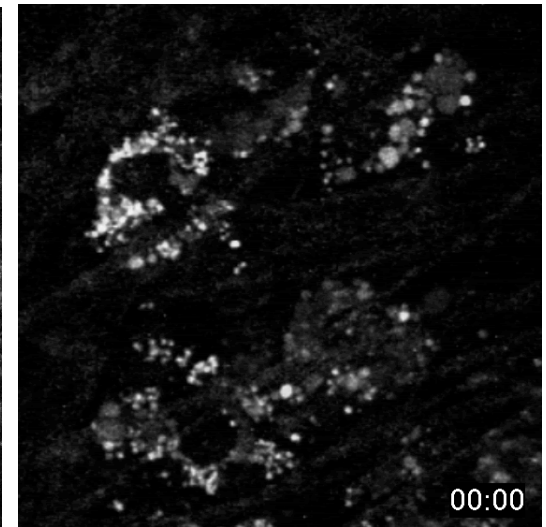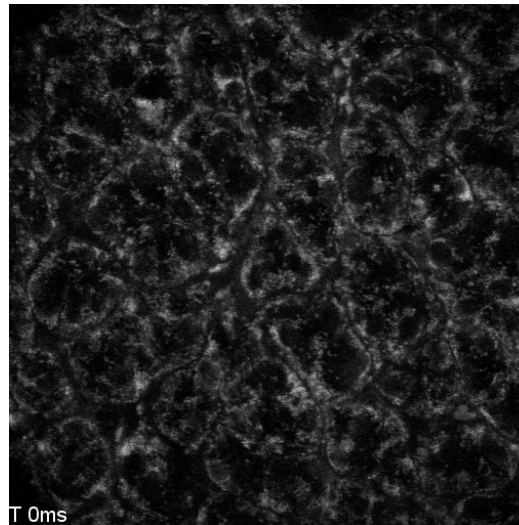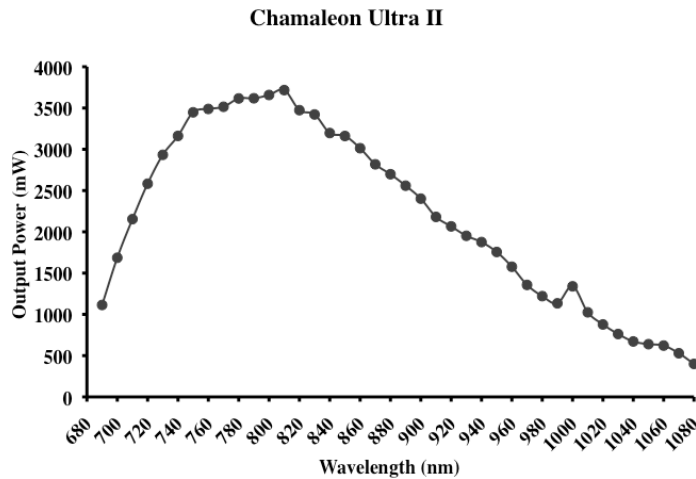

# Loss of power throughout the optics

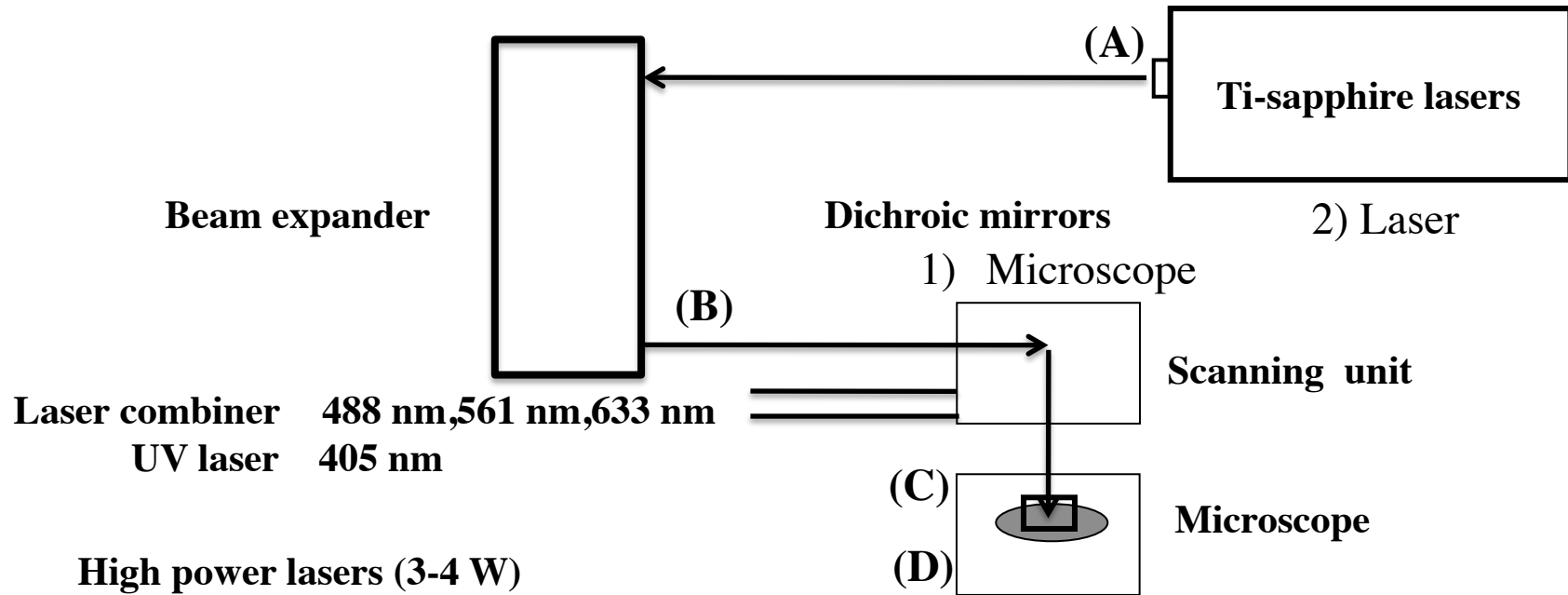


Chameleon Ultra II

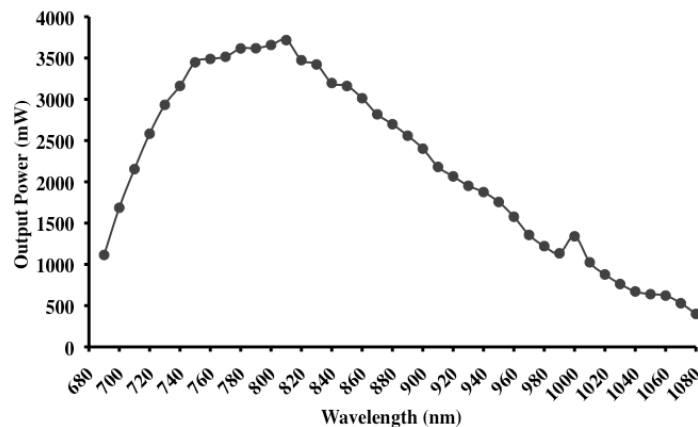


- (A) 800 nm – 3520 mW (100%)
- (B) 3280 mW (93%)
- (C) 864 mW (25%) – S.H. cuts 20%
- (D) 340 mW (10%) – 60x N.A 1.2


# Control the power at the specimen

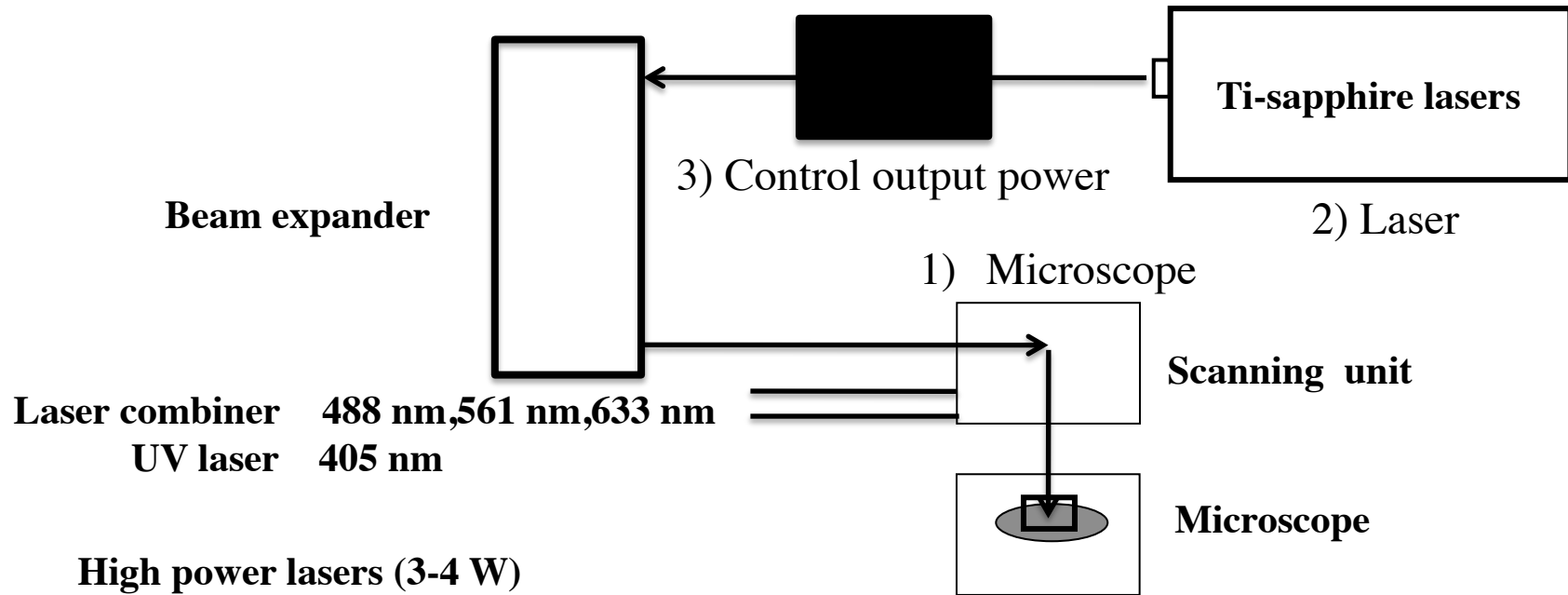

- 1) ND filters
  - a) Single
  - b) Carousel with multiple filters (8-10)
    - a) Manually or software controlled
- 2) ND continuous filter wheel

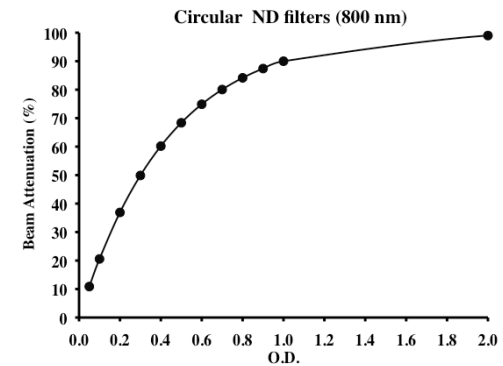

# Control the power at the specimen

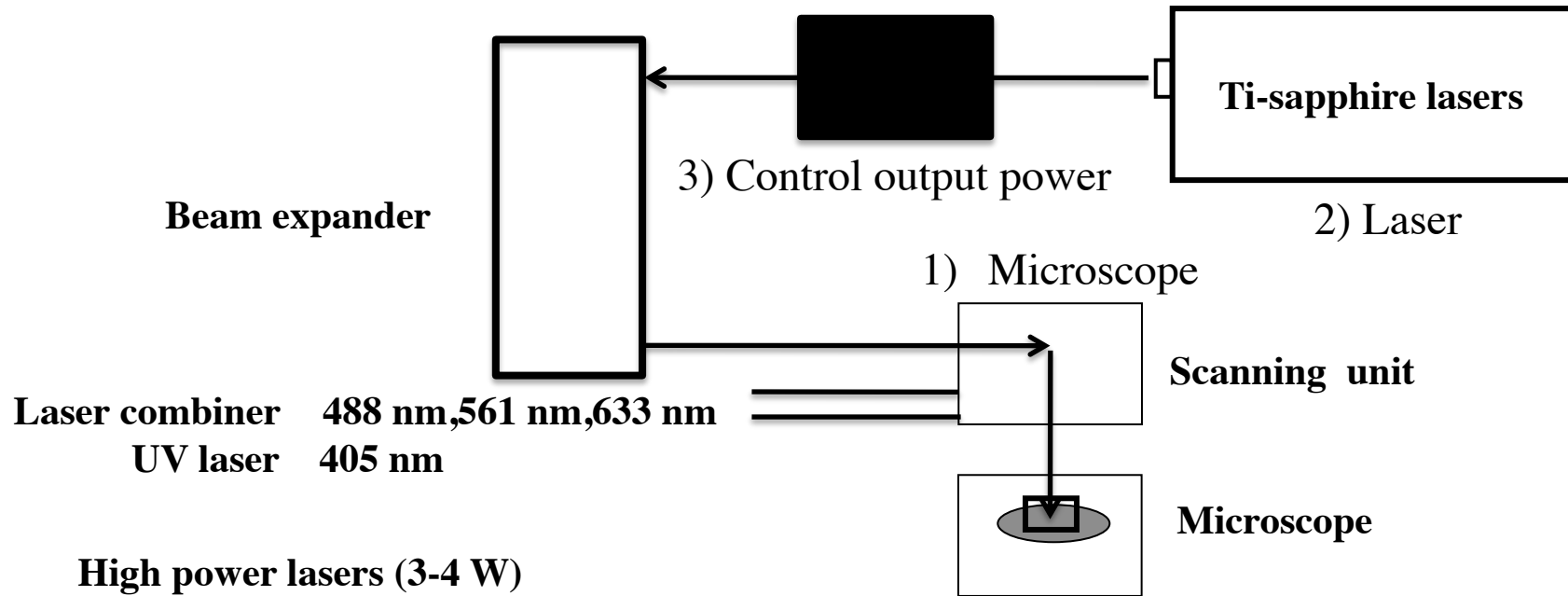

1) ND filters

2) ND continuous filter wheel

3) AOM (Acousto-optic modulator)

4) EOM (Electro-optic modulator)

a) Easy integration with the software

b) Size of the beam matching the aperture of the AOM

c) Significant pulse broadening (up to 600 fs)

a) Need for a pre-chirping system

d) Deflection of the beam

a) Not practical if different wavelengths are needed

b) Need for an automatic realignment set up (expensive)

# Broadening of the pulse width

(a)




Beam expander

Laser combiner 488 nm, 561 nm, 633 nm  
 UV laser 405 nm


Pulses: 100-150 fs

(a) Measure the pulse at the “source”


(b) Measure the pulse and the power at the “specimen”

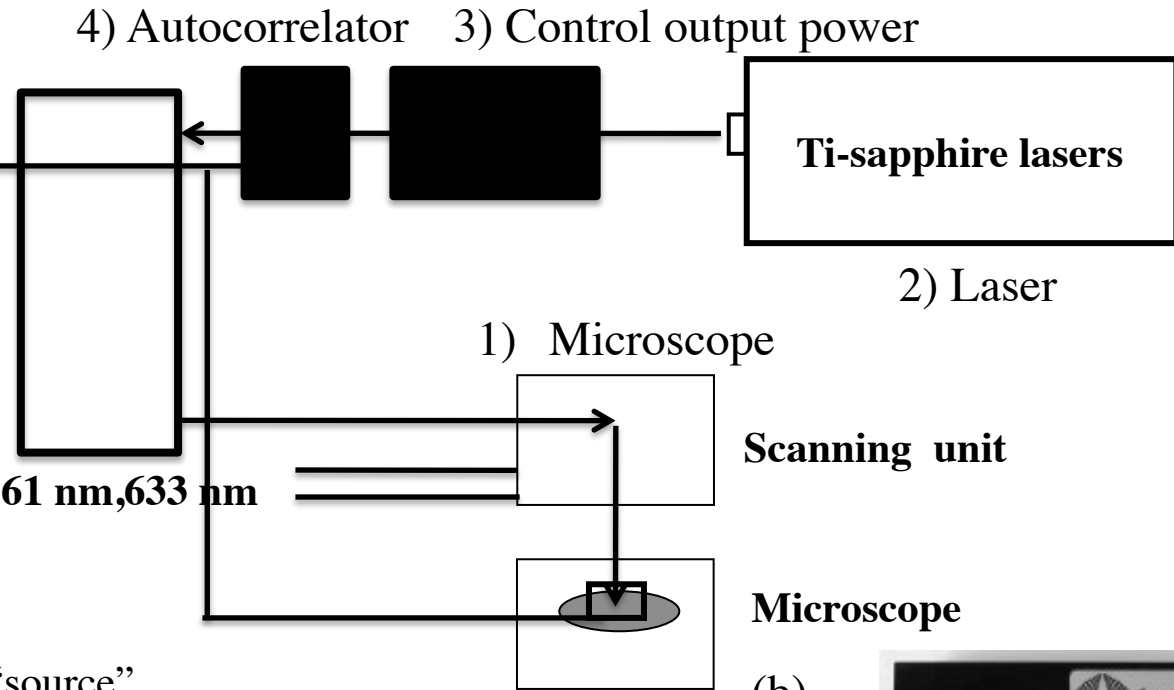

(b)

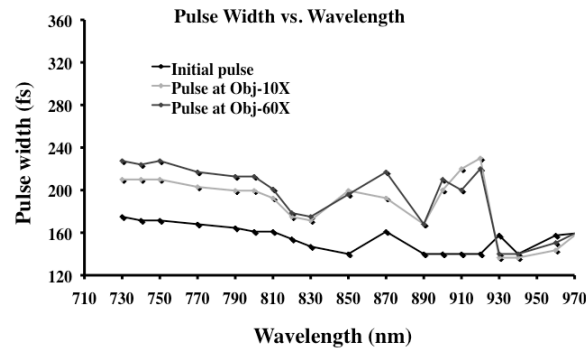

# Size of the laser beam

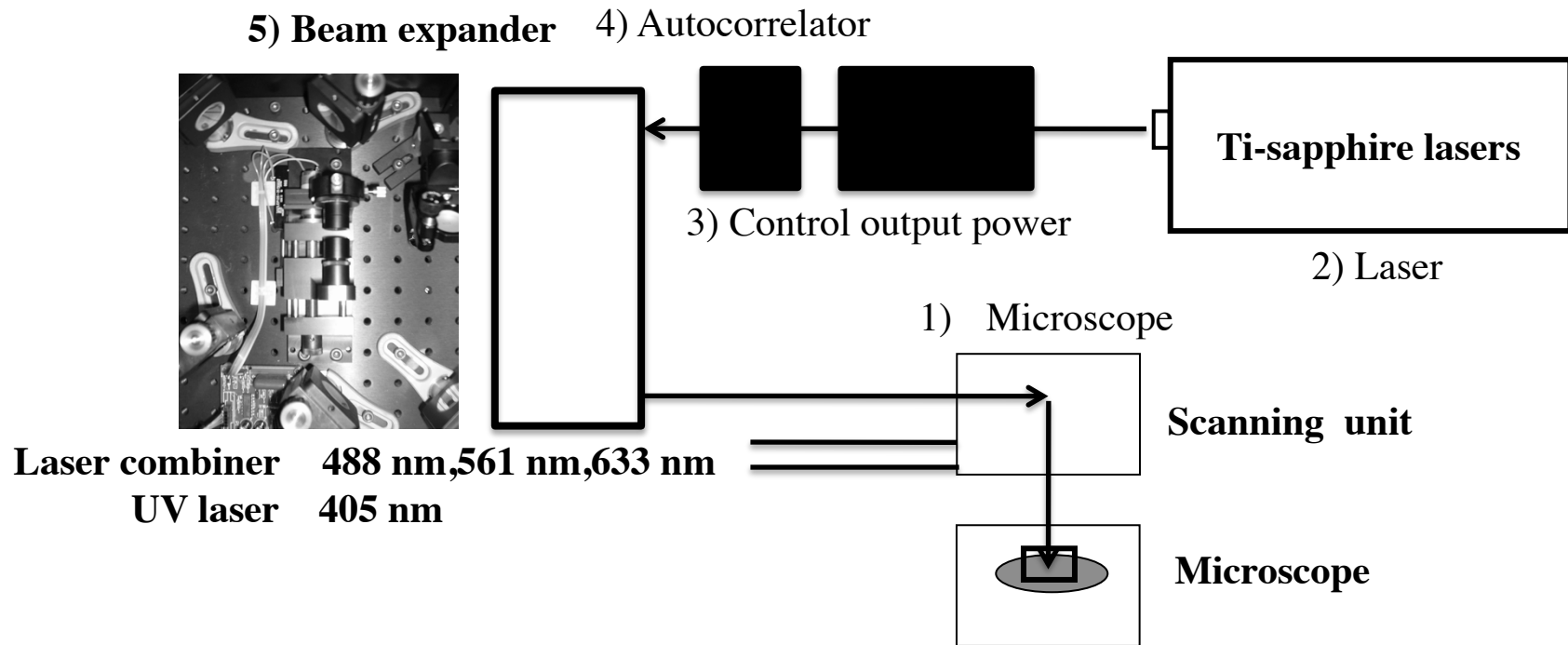


1) Control the size of the beam

2) Control the power at the specimen

# Filling the backaperture of the lens

## Filling the back aperture of the lens

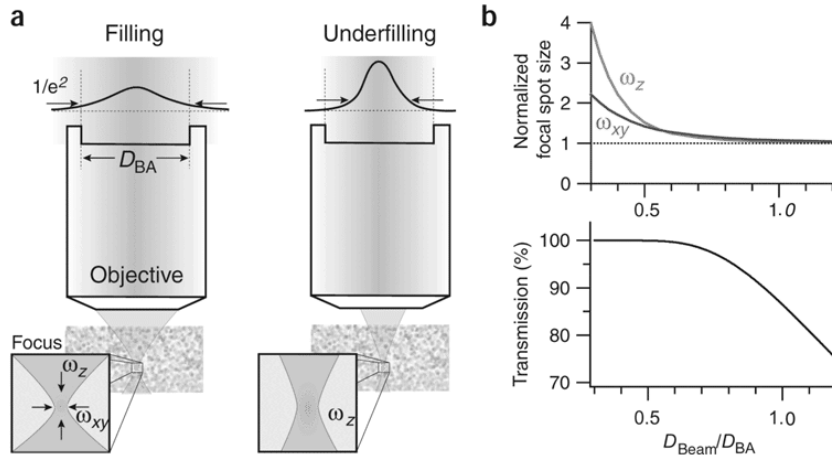


Essential for large lenses such as the 20X



## Underfilling





## Filling

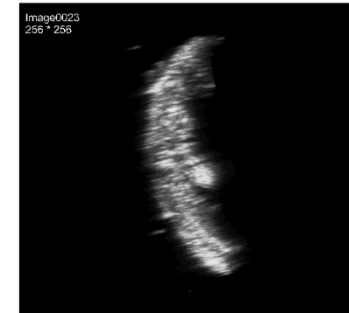

## Control the power at the specimen by overfilling

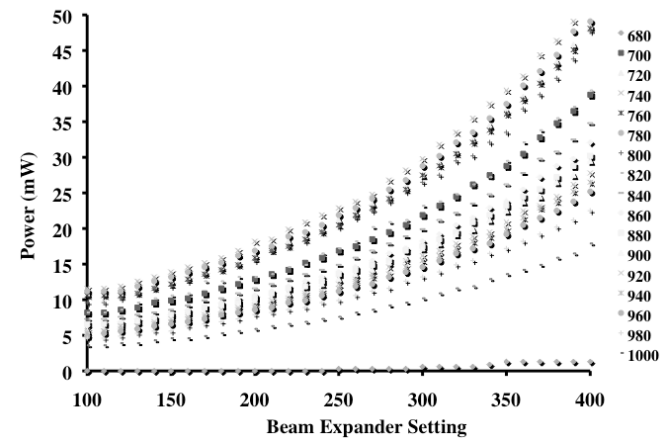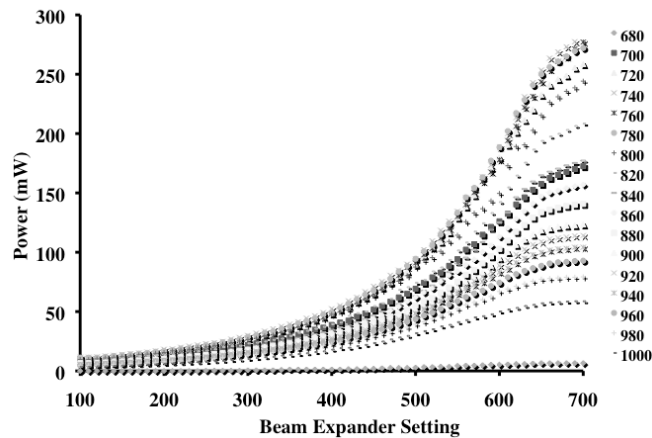


# Challenge: alignment of the beam

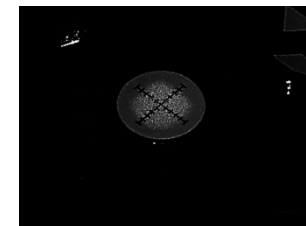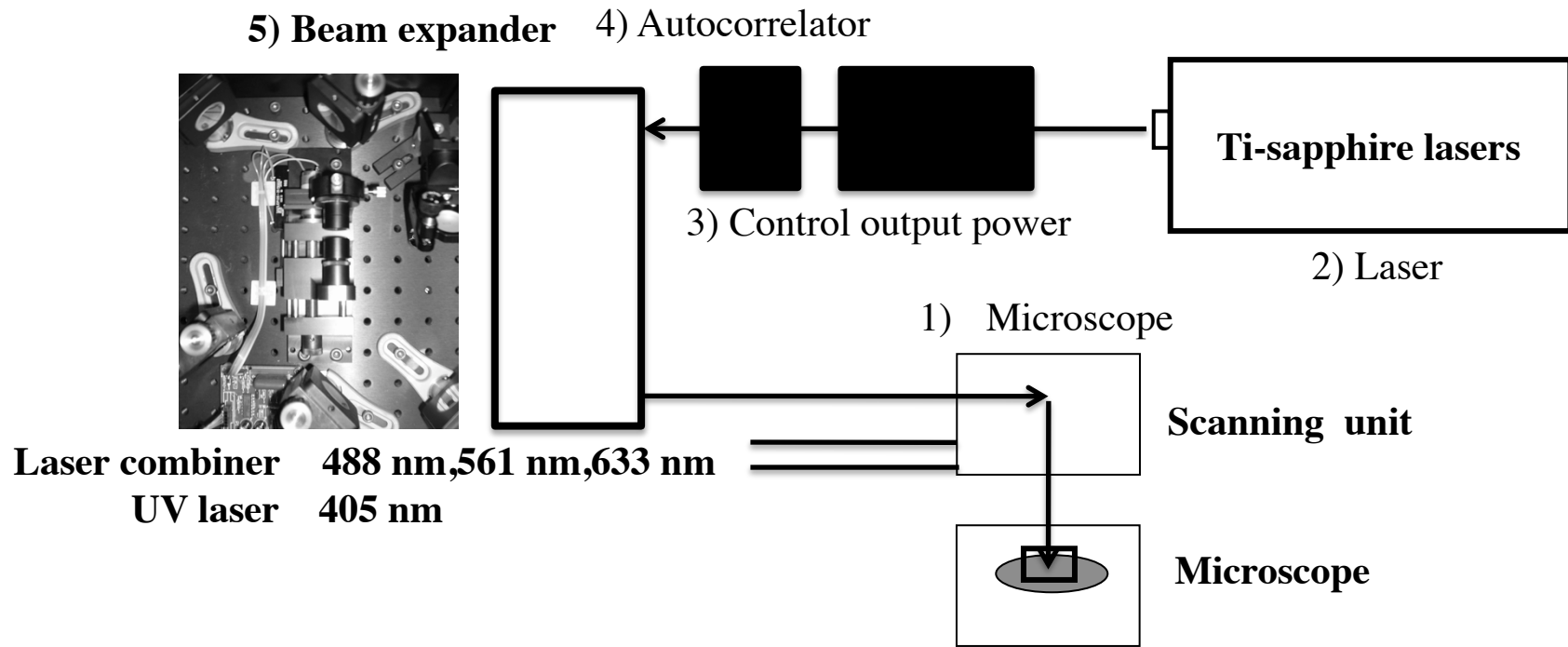




# Proper optics

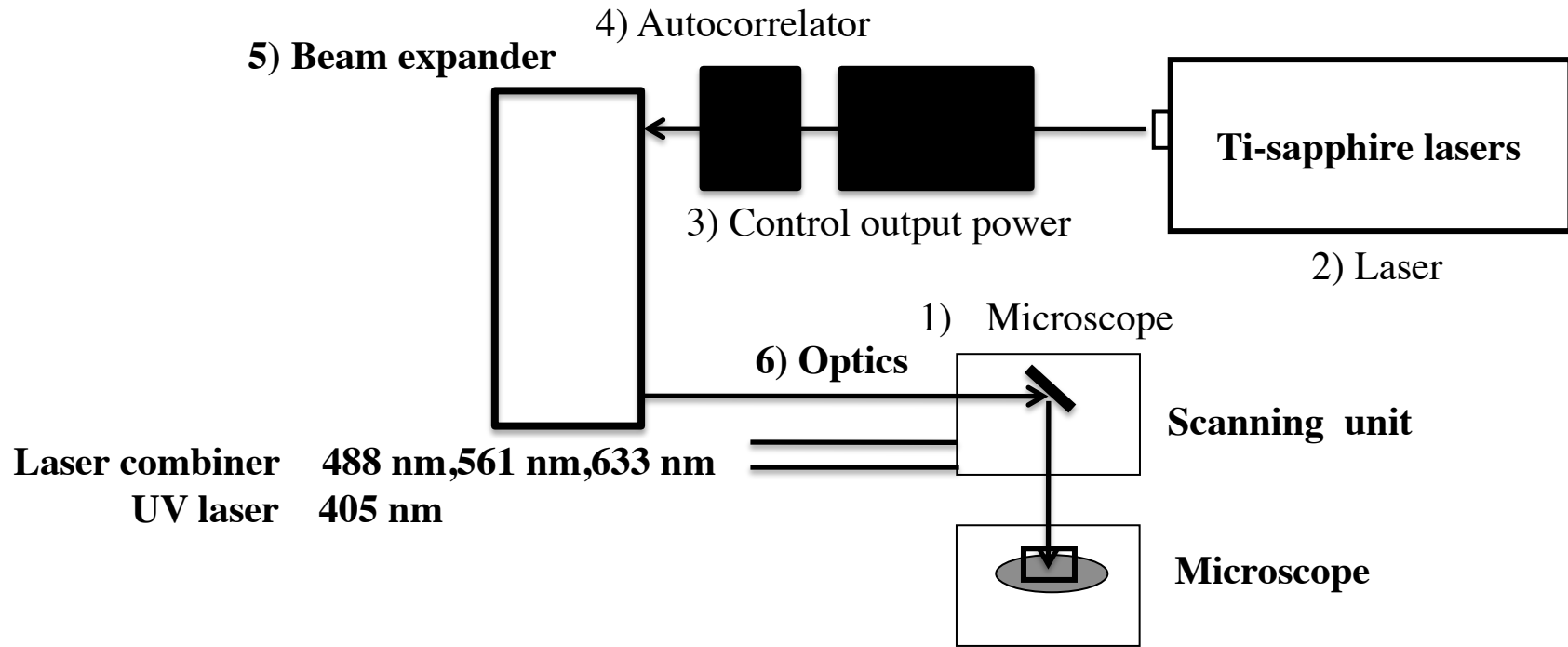


**6) Excitation Dichroic mirror – reflect above 675-680 nm**

# Non-descanned detectors

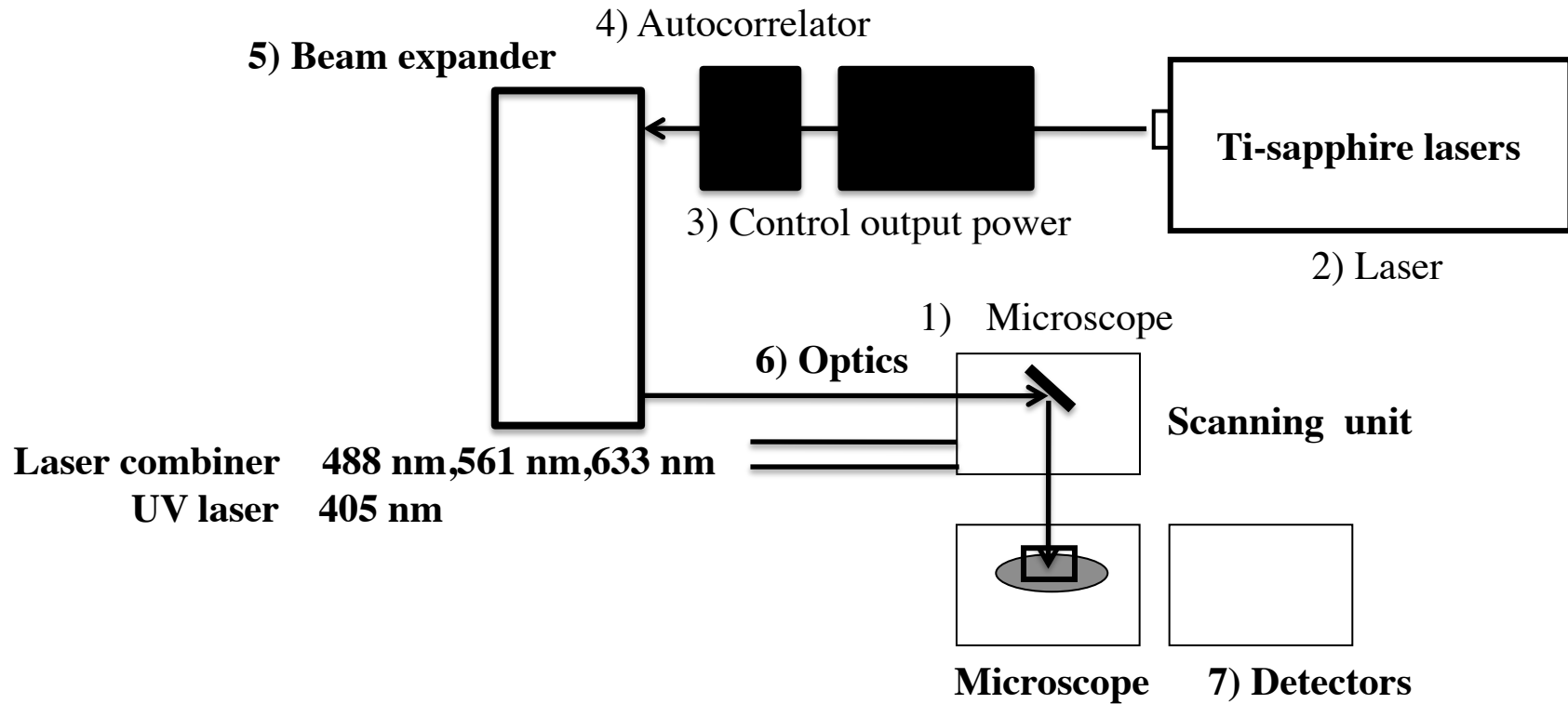



**A) Descanned detectors**

**B) Non-descanned detector**

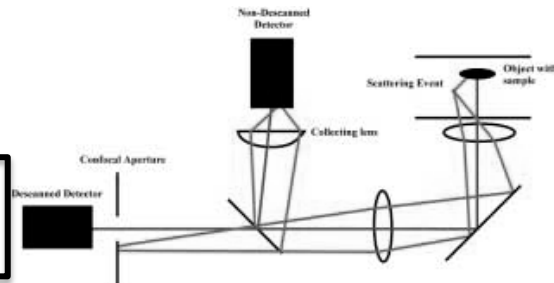


**Positioning**

# Non-descanned detectors

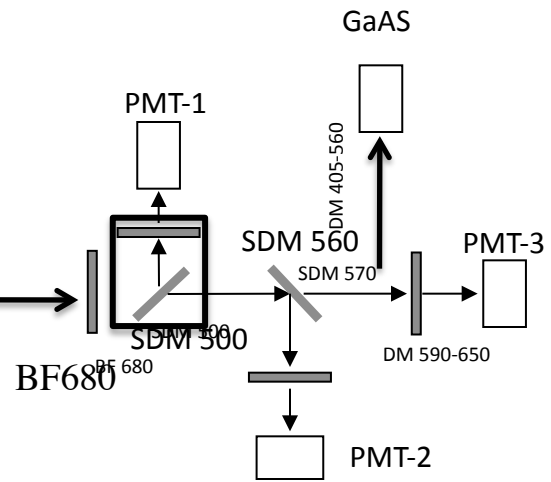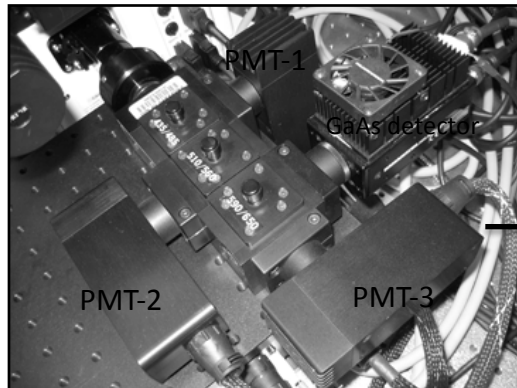

**3 Cooled PMT from Hamamatsu R6060-11**

**1 Gallium Arsenide PMT**

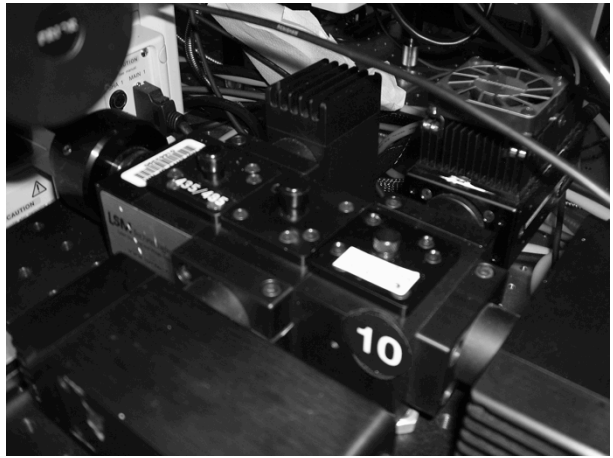


20x  
excitation 750  
emission < 510

Cooled PMT

GaAs detector

GaAs det gradient


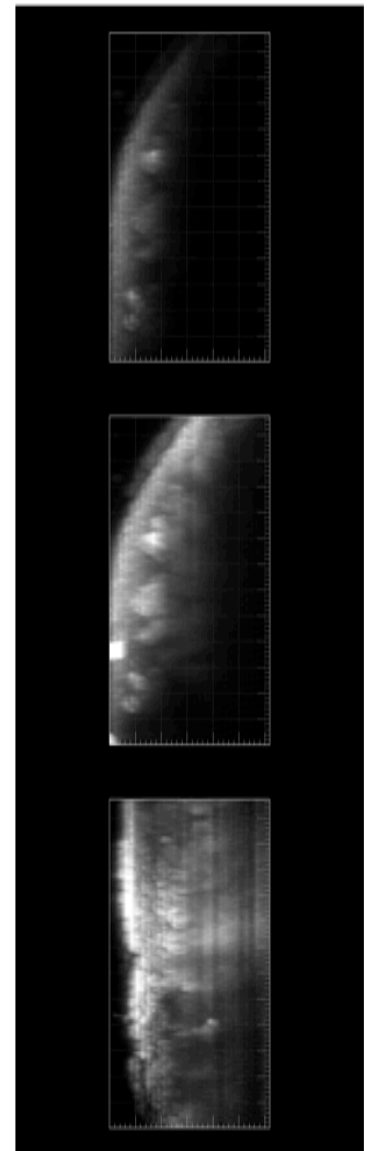


0 300  $\mu\text{m}$

# Non-descanned detectors

## Objective inverter with PMT



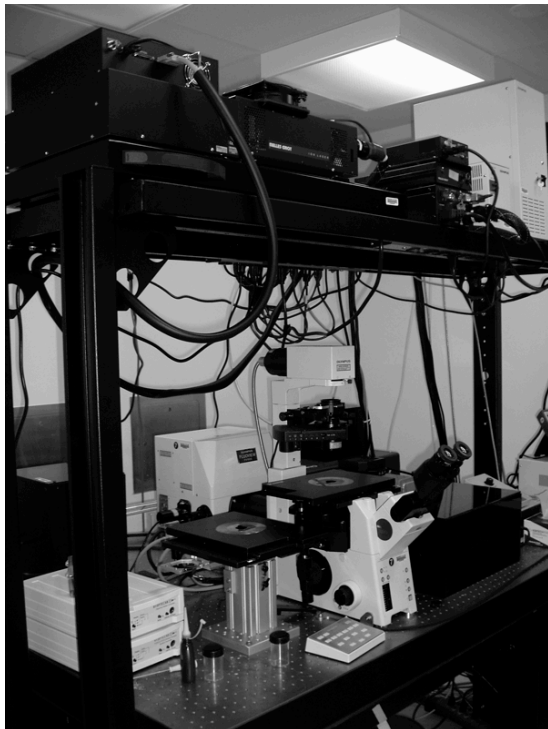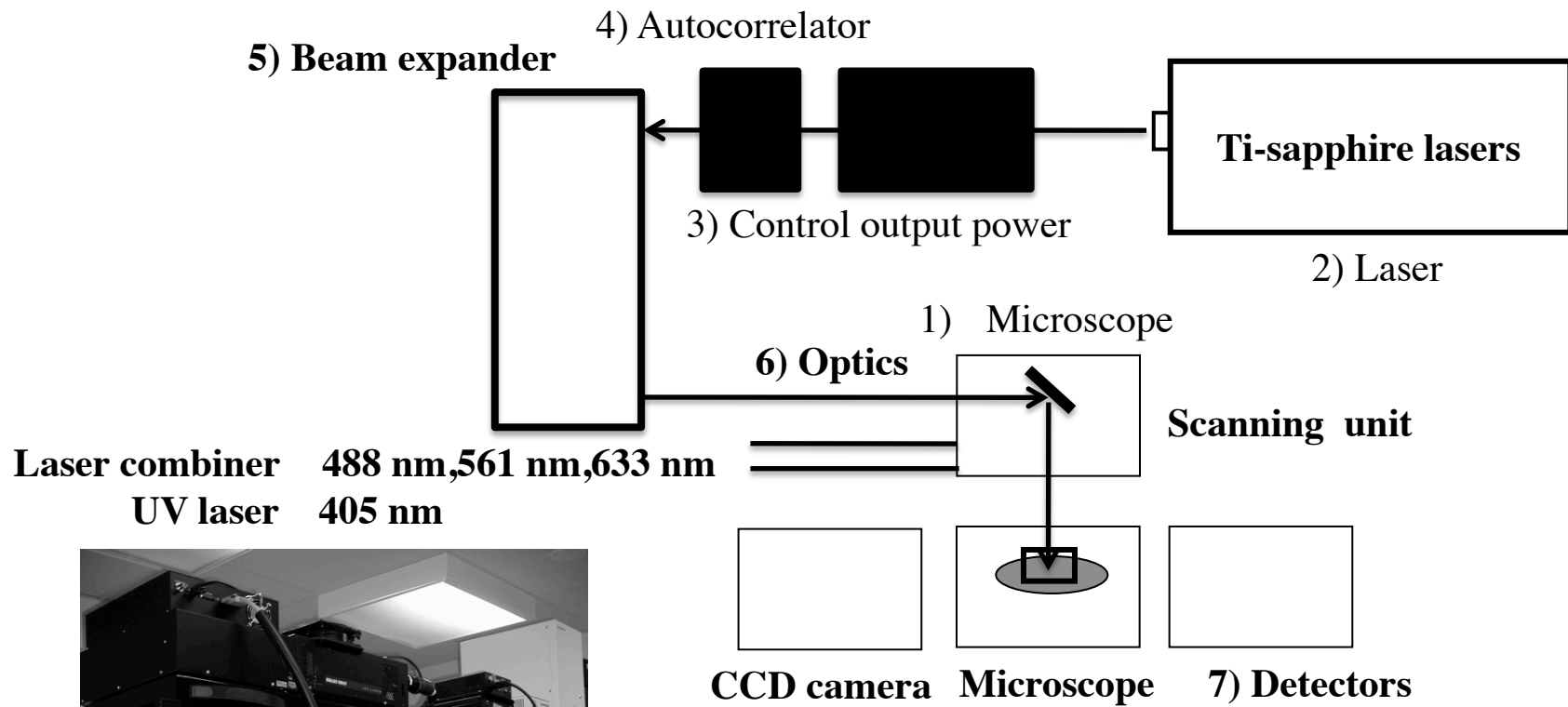

Confocal microscopy

Two-photon microscopy